



Visit the webpage below to
see supported devices/
download app (scan the QR
code for easy access).

www.i-sens.com/smartlog

i·sens



i-SENS, Inc.

43, Banpo-daero 28-gil
Seocho-gu, Seoul 06646, Korea
www.i-sens.com

Australian Distributor:

Pharmaco (Australia) Ltd
Chatswood, NSW 2067
Phone: 1800 114 610
www.pharmacodiabetes.com.au

NZ Distributor:

Pharmaco (NZ) Ltd
Auckland
Phone: 0800 GLUCOSE (0800 458267)
www.pharmacodiabetes.co.nz

© 2023 i-SENS, Inc. All Rights Reserved. PGUAA-0000292 REV1 2023-10   0123 

i·sens

User Manual

CareSens™ **N** **Premier** Blood Glucose Monitoring System

- Easy Data Transfer by Bluetooth Connection
- 2.5" Large, Backlit Display
- Test Result Tagging Options



Welcome to the CareSens N Premier Blood Glucose Monitoring System

Thank you for choosing the CareSen N Premier Blood Glucose Monitoring System. The system provides you with rapid and convenient blood glucose *in vitro* (i.e., outside the body) diagnostic monitoring. You can obtain accurate results in just 5 seconds with a small (0.5 µL) blood sample. Test results can be sent to a smartphone through Bluetooth wireless connection.

Bluetooth wireless technology

Bluetooth wireless technology is used by some smartphones and many other devices. CareSens N Premier Meter uses Bluetooth wireless technology to pair* and to send your glucose results to a smartphone.

The CareSens N Premier Meter is designed to work with the SmartLog mobile app. When using the CareSens N Premier System, we suggest you pair your CareSens N Premier Meter with your smartphone and track your results.

*The process of creating a connection between two Bluetooth devices requires an auto-generated passkey to be exchanged between the two devices. Once the devices are paired, they will automatically communicate with each other when the Bluetooth feature is activated.

- No part of this document may be reproduced in any form or by any means without the prior written consent of i-SENS.
- The information in this manual is correct at the time of printing. However, i-SENS reserves the right to make any necessary changes at any time without notice as our policy is one of continuous improvement.

 **WARNING:**

Turn off the Bluetooth feature in areas where the use of wireless devices is restricted, such as hospitals, some healthcare professional offices, and aeroplanes.

Trademarks

The Bluetooth® word mark and logos are registered trademarks owned by Bluetooth SIG, Inc. and any use of such marks by i-SENS, Inc. is under license. All other trademarks and trade names are those of their respective owners.

Note:

The SmartLog™ mobile app may not be compatible with all smartphones.



Table of Contents

Information

| | |
|--|----|
| Important Information: Read This First! | 6 |
| Specifications | 8 |
| CareSens N Premier Blood Glucose Monitoring System | 9 |
| Inserting or Replacing the Batteries | 10 |
| Caring for Your System | 13 |
| CareSens N Blood Glucose Test Strip | 14 |
| CareSens N Premier Blood Glucose Meter | 16 |
| CareSens N Premier Blood Glucose Meter Display | 17 |

Preparation

| | |
|--|----|
| Setting Up Your System | 18 |
| Setting Up Bluetooth | 19 |
| Adjusting the Date and Time | 22 |
| Setting the Sound On/OFF | 24 |
| Turning on the Strip Expiration Date Indicator | 25 |
| Turning on the Hypoglycemia (HYPo) Indicator | 26 |
| Setting the Strip Expiration Date Indicator | 27 |
| Checking the System | 28 |
| Control Solution Testing | 29 |
| Comparing the Control Solution Test Results | 32 |

Testing

| | |
|---------------------------------------|----|
| Using the Lancing Device | 33 |
| Preparing the Lancing Device | 34 |
| Preparing the Meter and Test Strip | 36 |
| Applying Blood Sample | 36 |
| Discarding Used Lancets | 40 |
| Alternate Site Testing | 41 |
| HI and Lo Messages | 44 |
| Target Blood Glucose Ranges | 45 |
| Transferring Test Results Using Cable | 45 |

Additional Functions

| | |
|---|----|
| Meter Memory | 46 |
| Viewing Averages Stored in Memory | 46 |
| Viewing Test Results Stored in Memory | 49 |
| Setting the Alarm Function | 50 |
| Setting the Post-meal Alarm (PP2 alarm) | 50 |
| Setting the Time Alarms (alarm 1–3) | 52 |

Maintenance

| | |
|---|----|
| Understanding Error Messages | 54 |
| General Troubleshooting | 56 |
| Performance Characteristics | 57 |
| System Accuracy and Measurement Precision | 57 |
| Influence Quantities | 58 |
| User Performance Evaluation | 60 |
| Warranty Information | 61 |

Important Information: Read This First!














For optimum benefits, please read the entire manual contents before using the system.

Intended use:

CareSens N Premier Blood Glucose Monitoring System is intended for the quantitative measurement of glucose in fresh capillary whole blood from the fingertip. Alternate site testing from the forearm and palm may also be used to measure glucose levels.

The CareSens N Premier Blood Glucose Monitoring System is intended for use outside the body (*in vitro* diagnostic use) and is intended for use as an aid to monitor the effectiveness of diabetes management. The system is for self-testing or healthcare professional use and should not be used for the diagnosis of or screening for diabetes or for testing newborns.

Meaning of Symbols Used:

-  For *in vitro* diagnostic use  Batch code
-  This product fulfills the requirements for Directive 98/79/EC on *in vitro* diagnostic medical devices
-  Caution, Warning (consult instructions for use for cautions and warnings)
-  Warning (Contains coin battery)
-  Use-by date
-  Do not discard this product with other household-type waste
-  Do not reuse  Manufacturer
-  Consult instructions for use  Serial number
-  Temperature limitation  EMC & Safety RCM Mark

- Glucose in blood samples reacts with the chemical in the test strip to produce a small electrical current. The CareSens N Premier meter detects this electrical current and measures the amount of glucose in the blood sample.
- The CareSens N Premier blood glucose meter is designed to minimise code related errors in monitoring by using the no-coding function.
- The CareSens N Premier blood glucose meter should be used only with the CareSens N Test Strips.
- An abnormally high or low red blood cell count (hematocrit level over 65 % or below 15 %) may produce inaccurate results.
- If your test result is below 3.3 mmol/L or above 13.3 mmol/L, consult a healthcare professional immediately.
- Inaccurate results may occur in severely hypotensive individuals or patients in shock. Inaccurate low results may occur for individuals experiencing a hyperglycemic-hyperosmolar state, with or without ketosis. Critically ill patients should not be tested with blood glucose meters.
- Inaccurate results may occur in patients undergoing oxygen therapy.

If you need assistance, please contact your authorised i-SENS sales representative or visit www.i-sens.com for more information.

Specifications

Product specifications

| | |
|-----------------------|--|
| Measurement range | 1.1–33.3 mmol/L |
| Sample size | Minimum 0.5 µL |
| Test time | 5 seconds |
| Sample type | Fresh capillary whole blood |
| Calibration | Plasma-equivalent |
| Assay method | Electrochemical |
| Battery life | 1,000 tests |
| Power | Two 3.0 V lithium batteries (disposable, type CR2032) |
| Memory | 1,000 test results |
| Size | 103 x 54 x 17 mm |
| Weight | 71.7 g (with batteries) |
| Bluetooth® technology | <ul style="list-style-type: none">• Frequency range: 2.4–2.4835 GHz• Operating range distance: maximum 10 meters (unobstructed)• Operating channels: 40 channels• Security encryption: 128-bit AES (Advanced Encryption Standard) |

Operating ranges

| | |
|-------------------|---------------------|
| Temperature | 5–50 °C (41–122 °F) |
| Relative humidity | 10–90 % |
| Hematocrit | 15–65 % |

Storage/Transport conditions

| | | |
|-------------------|--------------------------------|---------|
| Temperature | Glucose meter (with batteries) | 0–50 °C |
| | Test strip | 1–30 °C |
| | Control solution | 8–30 °C |
| Relative humidity | Test strip | 10–90 % |

CareSens N Premier Blood Glucose Monitoring System

CareSens N Premier Blood Glucose Monitoring System includes the following items:

- * CareSens N Premier Blood Glucose Meter
- * User Manual
- * Quick Reference Guide
- * Batteries
- * CareSens N Blood Glucose Test Strips
- * Lancets
- * Lancing Device
- * Logbook
- * Carrying Case
- * Registration card

- Check all the components after opening the CareSens N Premier Blood Glucose Monitoring System package. The exact contents are listed on the main box.
- The cable for data management software can be ordered separately. Please contact your authorised i-SENS sales representative.

Inserting or Replacing the Batteries



WARNING

- Chemical burn hazard.
- Keep batteries away from children.
- This product contains lithium button/coin cell batteries. If a new or used lithium button/coin cell battery is swallowed or enters the body, it can cause severe internal burns and can lead to death in as little as 2 hours. Always completely secure the battery compartment. If you think batteries might have been swallowed or placed inside any part of the body, seek immediate medical attention.
- The cells shall be disposed of properly, including keeping them away from children, and even used cells may cause injury.
- If you suspect your child has swallowed a button/coin cell battery immediately call the Poisons Information Centre on 13 11 26 in Australia, or 0800 764 766 in New Zealand.

The CareSens N Premier meter uses two CR2032 lithium batteries. Before using the meter, check the battery compartment and insert batteries if empty. When the **+** **-** symbol appears on the display while the meter is in use, the batteries should be replaced as soon as possible. The test results may not be saved if the batteries run out.

Step 1





Make sure the meter is turned off.
To open the battery cover, insert a narrow object, such as a ballpoint pen into the slot, and lift the cover up.

Note:

Please take care when opening the battery compartment – the cover may spring up.

Step 2

Remove the used batteries one at a time. Slip your index finger under the battery to lift and pull out as shown. Insert two new batteries with the + side facing up and make sure the batteries are inserted firmly.

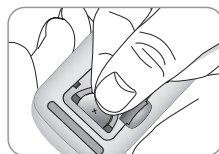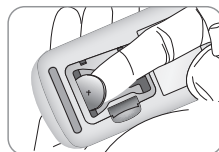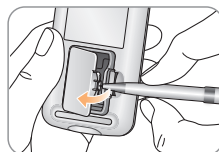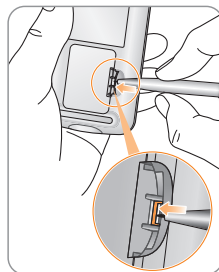

Step 3

Place the cover on the battery compartment, and push it down until it clicks into place.



Note:

Removing the meter batteries will not affect your stored results. However you may need to reset your meter settings. See page 18.

Caring for Your System

Use a soft cloth or tissue to wipe the meter exterior. If necessary, dip the soft cloth or tissue in a small amount of alcohol. Do not use organic solvents such as benzene or acetone, or household and industrial cleaners that may cause irreparable damage to the meter.

Caution:


- Do not expose the meter to direct sunlight, heat, or excessive humidity for an extended period of time.
- Do not let dirt, dust, blood, or water enter into the meter's test strip port.
- Avoid getting any liquid or moisture in the test strip vial. This can affect the test strips and cause inaccurate test results.
- Do not apply samples other than capillary whole blood or control solution to the test strip.
- Do not drop the meter or submit it to strong shock.
- Do not try to fix or alter the meter in any way.
- Strong electromagnetic radiation may interfere with the proper operation of this device. Keep the device away from sources of strong electromagnetic radiation, especially when measuring your blood glucose.
- Store all the meter components in the carrying case to prevent loss and help keep the meter clean.

Disposal of the meter

If you need to throw your meter away, you should follow existing policies and procedures of your own country or region. For information about correct disposal, please contact your local council or authority. If you need assistance, contact your authorised i-SENS sales representative or visit www.i-sens.com.

CareSens N Blood Glucose Test Strip

The CareSens N Premier Blood Glucose Monitoring System measures blood glucose quickly and accurately. It automatically absorbs the small blood sample applied to the narrow edge of the strip.

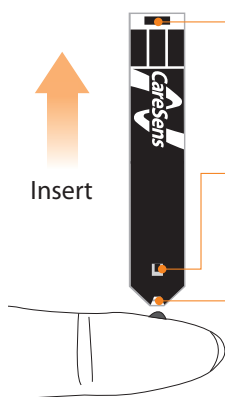

Contact bars

Gently push the test strip, with its contact bars facing up, into the test strip port of meter

Confirmation window

Check here to see whether sufficient blood sample has been applied

Edge to apply blood sample

Apply blood sample here for testing

Caution:

- The CareSens N Test Strips should be used only with fresh capillary whole blood samples.
- Do not reuse test strips.
- Do not use test strips past the expiration date.
- Test strips in new, unopened vials and test strips in vials that have been opened can be used up until the expiration date printed on the test strip box and vial label if the test strips are used and stored according to its storage and handling methods.
- Store test strips in a cool and dry place at a temperature between 1–30 °C (34–86 °F).
- Keep test strips away from direct sunlight or heat and do not freeze.
- Store test strips only in their original vial.
- Close the vial tightly after taking out a test strip for testing and use the strip immediately.
- Handle test strips only with clean and dry hands.
- Do not bend, cut, or alter test strips in any way.
- For detailed storage and usage information, refer to the CareSens N test strip package insert.
- Keep the meter and testing supplies away from young children.
- Drying agents in the vial cap may be harmful if inhaled or swallowed and may cause skin or eye irritation.

CareSens N Premier Blood Glucose Meter

Data Port

Used to transfer data from the meter to a computer with a cable

Display

Shows results, messages

◀, ▶ Button

Turns the meter on, selects or changes information

S Button

Turns the meter on/off, confirms menu selections, and changes information

Test Strip Port

Insert test strip here

Test Strip Ejector

Slide down to discard the used strips



Note:

- The cable for data management software can be ordered separately. Please contact your authorised i-SENS sales representative.
- The unit of measurement is fixed to mmol/L and it cannot be changed by the user.

CareSens N Premier Blood Glucose Meter Display



- ① **Test results:** test results displaying panel
- ② **Memory Symbol:** appears when test results stored in the memory are displayed
- ③ **PP2 alarm:** appears when the post-meal alarm has been set
- ④ **Bluetooth symbol**
- ⑤ **Mute symbol:** appears only when the sound is set to OFF
- ⑥ **Control Solution flag:** appears when the control solution test results are saved or displayed

- ⑦ **Decimal point:** appears when the blood glucose measuring unit is set to mmol/L
- ⑧ **mmol/L, mg/dL:** unit for measuring blood glucose
- ⑨ **alarm:** appears when the time alarm has been set
- ⑩ **Battery symbol:** indicates meter battery is running low and needs to be replaced
- ⑪ **Blood insertion symbol:** indicates meter is ready for the application of a drop of blood or control solution
- ⑫ **Pre-meal test flag:** used for tests done before eating
- ⑬ **Post-meal test flag:** used for tests done after eating
- ⑭ **Fasting test flag:** used for tests done after fasting for at least 8 hours
- ⑮ **Month/Day/Hour/Minute:** appears date and time

Note:


It is recommended to check if the display screen on the meter matches the illustration above every time the meter turns on. Do not use the meter if the display screen does not exactly match the illustration as the meter may show incorrect results.

Setting Up Your System

Press and hold the **S** button for 3 seconds to enter SET mode. After all settings are finished, press and hold the **S** button for 3 seconds to turn off the meter. Press the ◀ or ▶ button to change values. Press and hold the ◀ or ▶ button to scroll faster.

Step 1 Entering the SET Mode

Press and hold the **S** button for 3 seconds to enter SET mode. After all the segments flash across the screen, 'SET' will show up. Press the ◀ or ▶ button to select 'YES' and press the **S** button to go to the next step.

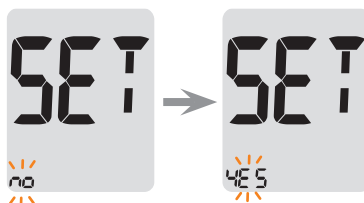

Note:

Follow steps 2–3 to pair your meter and smartphone. Pairing allows the meter to communicate wirelessly with your smartphone. Ensure that devices are within the maximum Bluetooth range (10 metres). Before pairing your meter and smartphone, download and install the SmartLog™ mobile app on your smartphone.

Setting Up Bluetooth

Step 2 Bluetooth Pairing


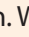
- ① If you do not want to connect your meter to your smartphone, press the **S** button when the screen shows on the right. The meter will go to Step 4 Year Setting mode.



- ② Press the ◀ or ▶ button. The meter screen will display 'OFF', 'On', and 'PAIR' in turn. To pair your meter with your smartphone, press the **S** button when 'PAIR' blinks on the screen.




Note:

The  symbol will appear on the screen when the Bluetooth feature is on. When the  symbol is not present on the screen, the Bluetooth feature is off. When you need to turn off/on the Bluetooth feature, press the **S** button when OFF/On blinks on the screen.



Step 3 Entering the PIN number

- 1 The  symbol and 'PIN' will appear if you press the **S** button when the meter screen shows 'PAIr'.




- 2 Launch SmartLog™ mobile app to start pairing the meter with your smartphone. Find the 'Accessories' menu and select the appropriate meter model on SmartLog™ mobile app.

Note: Follow the on-screen instructions on SmartLog™ mobile app to begin searching your meter. Some content or menus may differ depending on your smartphone's operating system or SmartLog™ version.

- 3 Look for 'CareSens' and the last 4 characters of the meter serial number on the SmartLog™ mobile app screen to correctly identify your meter. Touch your meter's ID (CareSens XXXX) on the SmartLog™ mobile app screen.

Serial
number

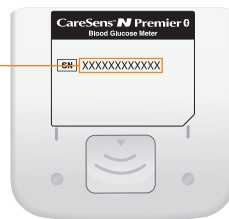

- 4 The meter will display six digit PIN number.



- ⑤ Enter the PIN number into the SmartLog™ mobile app and touch 'OK'. Make sure the PIN you enter on your smartphone matches the PIN on your meter screen.

- ⑥ When your meter and smartphone are paired and connected, the meter will display 'SUCCESS' and the saved test results will be transferred to your smartphone.



- ⑦ When the data transfer is finished, the meter will display 'End' on the screen. Press the **S** button to go to Step 10 Sound Setting mode. See page 24. If the meter displays 'FAIL' and then 'OFF', repeat steps 2 to 5.



Note:

Some smartphones, especially those that are not tested or approved by i-SENS, may be incompatible with your meter. Visit www.i-sens.com/smartlog for more information about supported smartphones. You can also scan the QR code on the back cover of this user manual.

Adjusting the Date and Time

Step 4 Setting the Year

Press the ◀ or ▶ button to adjust until the correct year appears. When the present year appears, press the **S** button to confirm your selection and to go to the next step.



Step 5 Setting the Month

A number indicating the month will blink on the screen.

Press the ◀ or ▶ button until the correct month appears. Press the **S** button to confirm your selection and to go to the next step.



Step 6 Setting the Date

Press the ◀ or ▶ button until the screen displays the correct date. Press the **S** button to confirm the date and to go to the next step.



Step 7 Setting the Time Format

The meter can be set in the AM/PM 12-hour or the 24-hour clock format.

Press the ◀ or ▶ button to select a format. The AM•PM symbol is not displayed in the 24-hour format. After selecting the format, press the **S** button to go to the next step.



Step 8 Setting the Hour

Press the ◀ or ▶ button until the correct hour appears. After the hour is set, press the **S** button to go to the next step.



Step 9 Setting the Minute

Press the ◀ or ▶ button until the correct minute appears. After setting the minute, press the **S** button to go to the next step.



Setting the Sound On/OFF

Step 10

On pressing the ◀ or ▶ button, the screen will display 'On' or 'OFF'. Press the **S** button to confirm the selection.


The meter will beep in the following instances if set to On.

- When you push a button to turn on the meter,
- When the test strip is inserted in the meter,
- When the blood sample is absorbed into the test strip and the test starts,
- When the test result is displayed,
- When you press and hold the ◀ button to set the post-meal (PP2) alarm,
- When it is time for a pre-set blood glucose test.



If the sound is set to OFF, none of the sound functions will work. After setting the sound, press the **S** button to progress to the next step.

Note:

The  symbol is displayed only when the sound is set to OFF.


Turning on the Strip Expiration Date Indicator

Step 11

This setting allows you to turn the strip expiration date indicator on or off. This setting turns the function on or off only. See page 27 to set the strip expiration date.


When 'EP' appears on the screen, press the ◀ or ▶ button. The screen will display 'On' or 'OFF'. Press the **S** button to confirm the setting.

If you do not want to set the indicator, press the **S** button while the screen displays 'OFF'.

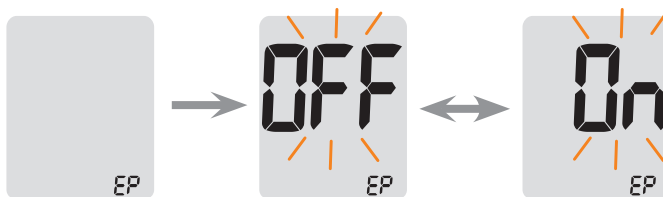

Note:


If the pre-set expiration date expires, the meter will display EP when the test strip is inserted. EP shows alternately also when the test result is displayed right after the test. If the expiration date is set to October of 2023, the meter will display EP at the start of November, 2023.

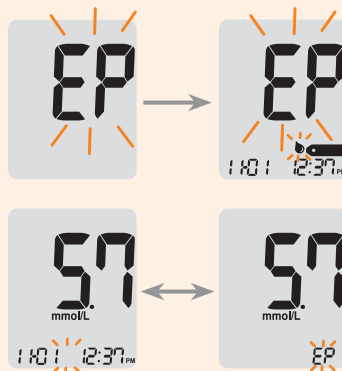

Turning on the Hypoglycemia (HYPo) Indicator


Step 12

This setting allows you to turn the hypoglycemia indicator (possible low blood sugar) On or OFF and to select the desired level for the indicator. You will be alerted any time your test result is lower than the selected level. On pressing the ◀ or ▶ button, the screen will display 'On' or 'OFF'. Press the **S** button when 'On' appears to enter the setting. Press the ◀ or ▶ button until the desired hypoglycemia level between 1.1 and 5.0 mmol/L appears. Then, press the **S** button to confirm the level and to return to step 2. See page 19.

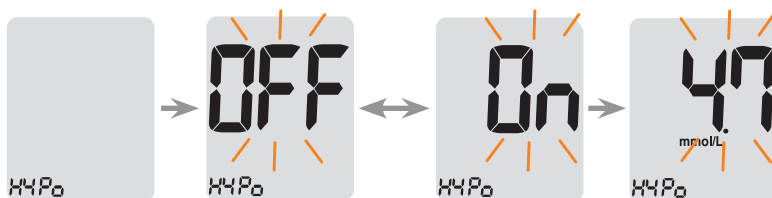

Note:

If the test result is lower than the pre-set hypoglycemia level, the meter will display the following.

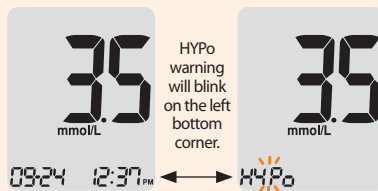

Caution:

Ask your healthcare professional to help you decide what your hypoglycemia level is before setting your level.

Setting the Strip Expiration Date Indicator

Step 1 Entering the Expiration Date Setting


Press and hold the ◀ and ▶ buttons at the same time for 3 seconds to enter the expiration date settings. After all segments flash across the screen, 'EP' will show up.

Note:

The strip expiration date is printed on the test strip vial.


Step 2 Setting the Year

A number indicating the year will blink in the left corner of the screen. Press the ◀ or ▶ button until the correct year appears. Press the **S** button to confirm the year and set the month.

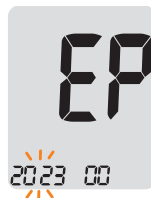

Step 3 Setting the Month

A number indicating the month will blink at the bottom of the screen. Press the ◀ or ▶ button until the correct month appears. After setting, press and hold the **S** button for 3 seconds to turn off the meter.

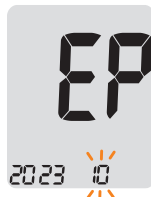

Checking the System



You may check your meter and test strips using the CareSens Control Solution (control A, B and/or C). The CareSens Control Solution contains a known amount of glucose and is used to check that the meter and the test strips are working properly. The test strip vials have CareSens Control Solution ranges printed on their labels. Compare the result displayed on the meter to the CareSens Control Solution range printed on the test strip vial. Before using a new meter or a new vial of test strips, you may conduct a control solution test following the procedure on pages 29-31.

Note:

- Use only the CareSens Control Solution (available for purchase separately).
- Check the expiration date printed on the bottle. When you first open a control solution bottle, record the discard date (date opened plus three (3) months) in the space provided on the label.
- Make sure your meter, test strips, and control solution are at room temperature before testing. Control solution tests must be done at room temperature (20–25 °C/68–77 °F).
- Before using the control solution, shake the bottle, discard the first few drops and wipe the tip clean.
- Close the control solution bottle tightly and store at a temperature between 8–30 °C (46–86 °F).

You may do a control solution test:

- When you want to practice the test procedure using the control solution instead of blood,
- When using the meter for the first time,
- Whenever you open a new vial of test strips,
- If the meter or test strips do not function properly,
- If your symptoms are inconsistent with the blood glucose test results and you feel that the meter or test strips are not working properly,
- If you drop or damage the meter.


Control Solution Testing

Step 1



Insert a test strip into the meter's test strip port, with the contact bars facing upwards.

Gently push the test strip into the port until the meter beeps. Be careful not to bend the strip while pushing it in.

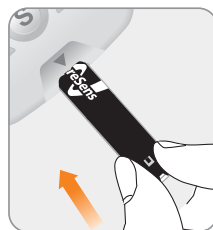The  symbol will show up.



Step 2


You can flag the control solution test result by pressing the  button for 3 seconds. To undo the control solution flag, press the  button for 3 seconds again.




Step 3


Shake the bottle before each test.

Remove the cap and squeeze the bottle to discard the first drop. Then wipe the tip with a clean tissue or cloth. Dispense a

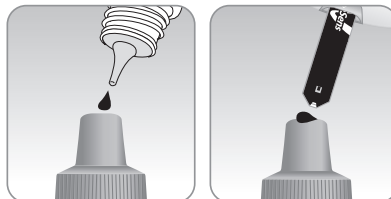

drop of control solution onto a clean non-absorbent surface. It helps to squeeze a drop onto the top of the cap as shown. After the  symbol appears on the display, apply the solution to the narrow edge of the test strip until the meter beeps. Make sure the confirmation window fills completely.


Note:

The meter may switch off if the control solution sample is not applied within 2 minutes of the  symbol appearing on the screen. If the meter turns off, remove the strip, reinsert, and start from step 1.

Step 4


The display segments will rotate clockwise and a test result will appear after the meter counts down from 5 to 1.

When flagged, the result is stored in the meter's memory but it is not included in the averages.

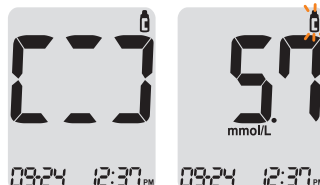

Step 5

Compare the result displayed on the meter to the range printed on the test strip vial. The result should fall within the range.

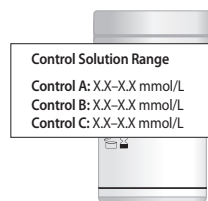

Caution:

The range printed on the test strip vial is for the CareSens Control Solution only. It has nothing to do with your blood glucose level.

Note:

The CareSens Control Solution can be ordered separately. Please contact your authorised i-SENS sales representative.

Comparing the Control Solution Test Results

The test result of each control solution should be within the range printed on the label of the test strip vial. Repeat the control solution test if the test result falls outside of the range. Out of range results may occur in following situations:


| Situations | Do This |
|---|---|
| <ul style="list-style-type: none">• When the control solution bottle was not shaken well,• When the meter, test strip, or the control solution were exposed to high or low temperatures,• When the first drop of the control solution was not discarded or the tip of the bottle was not wiped clean,• When the meter is not functioning properly. | Repeat the control solution test by referring to the Note on page 28. |
| <ul style="list-style-type: none">• When the control solution is past the expiration date printed on the bottle,• When the control solution is past its discard date (the date the bottle was opened plus three (3) months),• When the control solution is contaminated. | Discard the used control solution and repeat the test using a new bottle of control solution. |

If results continue to fall outside the range printed on the test strip vial, the CareSens N Test Strip and CareSens N Premier meter may not be working properly. Do not use your system and contact i-SENS sales representative.

Using the Lancing Device

You will need a lancing device in order to collect a blood sample.

You may use the lancing device that is included in the CareSens N Premier Blood Glucose Monitoring System or any other medically approved lancing device.

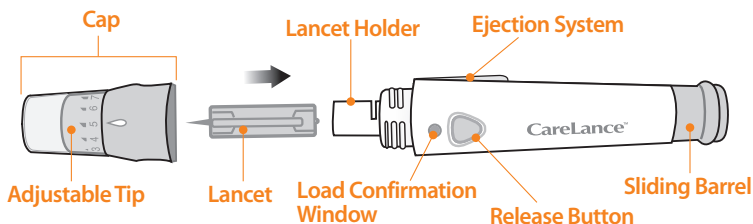

- The lancing device is for use by a single user only and should not be shared with anyone.
- Use a soft cloth or tissue to wipe the lancing device. If necessary, a small amount of alcohol on a soft cloth or tissue may be used.

Caution:

To avoid infection when drawing a sample, do not use a lancet more than once, and:

- Do not use a lancet that has been used by others.
- Always use a new sterile lancet.
- Keep the lancing device clean.

Note:

Repeated puncturing at the same sample site may cause pain or skin calluses (thick hard skin). Choose a different site each time you test.

Preparing the Lancing Device


Step 1

Wash hands and sample site with soap and warm water. Rinse and dry thoroughly.




Step 2

Unscrew and remove the lancing device tip.

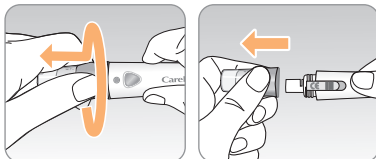


Step 3

Firmly insert a new lancet into the lancet holder. Hold the lancet firmly. Gently twist to pull off protective disk. Save disk to recap lancet after use. Replace lancing device tip.

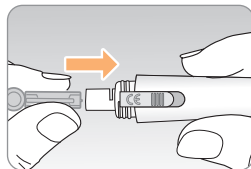

Step 4


Turn the adjustable tip until it is aligned with the load confirmation window and release button as shown.

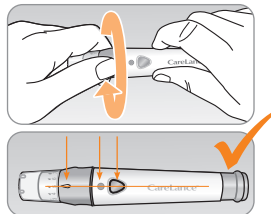

Step 5

The lancing device has seven puncture depth settings, numbered 1 through 7 (1 for a shallow puncture, 7 for a deeper puncture).

Choose a depth by rotating the top portion of the adjustable tip until the desired number aligns with the arrow.

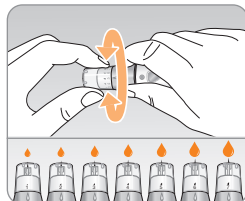


Note:

1 = least penetration of lancet into the skin.

7 = most penetration of lancet into the skin.

Step 6

To cock the lancing device, hold the body of lancing device in one hand and pull the sliding barrel with the other hand. The device is loaded when you feel a click and the load confirmation window turns red.

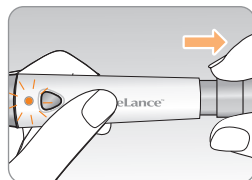


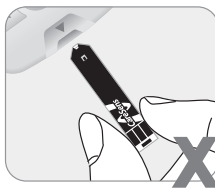
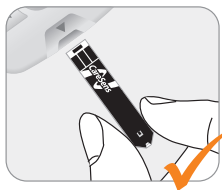
Note:

The skin depth to get blood samples will vary for various people at different sample sites. The lancing device's adjustable tip allows the best depth of skin penetration to get an adequate sample size.

Preparing the Meter and Test Strip


Step 7


Insert a test strip with the contact bars facing upwards into the meter's test strip port. Push the strip in gently until the meter beeps. Be careful not to bend the test strip. The  symbol will appear on the screen.




Applying Blood Sample

Step 8


Obtain a blood sample using the lancing device. Place the device against the pad of the finger. The best puncture sites are on the middle or ring fingers. Press the release button. Remove the device from the finger. Wait a few seconds for a blood drop to form. A minimum volume of 0.5 microliter is needed to fill the confirmation window (actual size of 0.5 μL : ).

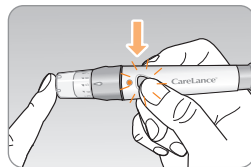


Step 9

After the  symbol appears on the screen, apply the blood sample to the narrow end of the test strip till the meter beeps. If the confirmation window is not filled in time because of abnormal viscosity (thickness and stickiness) or insufficient volume, the **Er4** message may appear.

It is recommended to place the test strip vertically into the blood sample site as shown below.

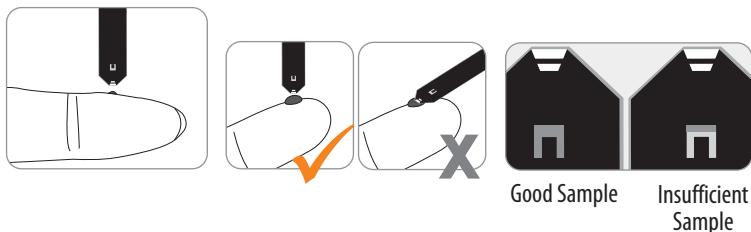


Note:

The meter may switch off if the blood sample is not applied within 2 minutes of the  symbol appearing on the screen. If the meter turns off, remove the strip and reinsert it, and start from Step 2.

Step 10

Apply the blood sample to the narrow end of the test strip until you hear a 'beep'. At this time, the display segments will rotate clockwise while the blood is going in.

The test result will appear after the meter counts down from 5 to 1. The result will be automatically stored in the meter's memory. If the test strip is removed after the test result is displayed, the meter will automatically switch off after 3 seconds. Discard used test strips safely in disposable containers. If the Bluetooth feature is activated, the meter will send the test result to the connected smartphone.

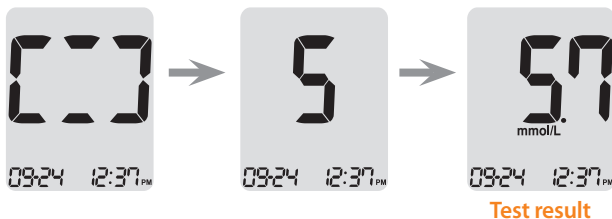

Note:

To transmit glucose data using the Bluetooth feature,

- The Bluetooth feature on the meter must be turned on,
- The meter and a smartphone must be paired,
- The SmartLog™ mobile app must be launched.

The meter will transmit data in the following cases,

- When the strip is ejected after measuring,
- When the meter is turned on (only when untransmitted data exists).




Step 11


You can attach a flag to a result to indicate particular situations while the strip is still in the meter. When the result is displayed right after a test, press the ◀ or ▶ button to select a pre-meal flag (🍏), a post-meal flag (🌞), a fasting flag (🕒), or a control solution flag (💧). When you remove the test strip while the desired flag is blinking, the test result is stored with the flag. If you do not want to add any flags on the test result, remove the strip after the test result is displayed.




No flag

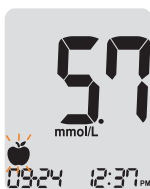

Pre-meal flag

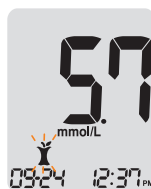

Post-meal flag

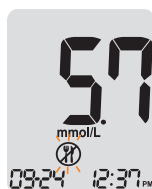

Fasting flag




Control solution flag

Discarding Used Lancets



Step 1

Unscrew the lancing device tip.

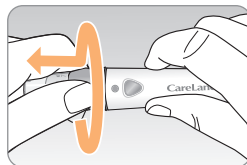

Step 2

Stick the lancet into the saved protective disk.
Push the lancet ejector forward with the thumb to dispose of the used lancet in a proper biohazard container.

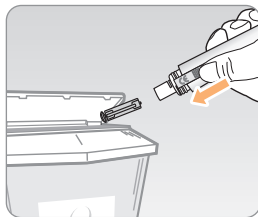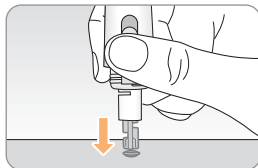

Caution:

The lancet is for single use only. Never share or reuse a lancet. Always dispose of lancets properly.


Alternate Site Testing

Before using AST, please consult your healthcare professional.

What is AST (Alternate Site Testing)?

Usually, we take the blood sample from the tip of the finger. However, since there are many nerve endings in the fingertip, it can be quite painful. When doing a glucose test, using different parts of the body such as the forearms and palms can reduce the pain during testing. This method of testing with different parts of the body is called Alternate Site Testing. While AST may reduce the pain during testing, it may not be simple for everyone and the following precautions should be observed during testing.

Alternate Sites for Testing

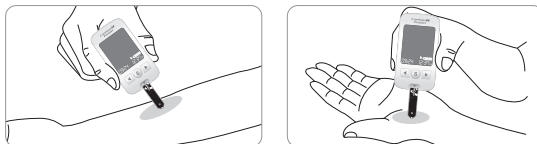

Alternate Site Blood Sampling (forearm and palm)

Select a clean, soft and fleshy sample site area free of visible veins and hair and away from bones. Gently massage the sample site to help blood circulation to minimise result differences between fingertip and alternate site sampling. Firmly press and hold the lancing device against site. Wait until the skin surface under the lancing device changes color. Then press the release button while continuing to apply pressure. Keep holding the lancing device against your skin until sufficient (at least 0.5 μ L, actual size: ●) blood is drawn. Carefully lift the lancing device away from your skin.

Things to Know When Using AST

Please read the following before testing at alternate sites (forearms and palms).

The capillary whole blood of the fingertips reflects changes in glucose levels more rapidly than in alternate sites. The test results from the fingertip testing and AST may differ due to factors such as lifestyle and ingested food which affect glucose levels.

Note:

Acceptable Situations for AST

When your blood glucose levels are stable

- Fasting period
- Before a meal
- Before sleeping

Situations Requiring Fingertip Test

When your blood glucose levels are unstable

- During the two hours after a meal or exercise
- When sick or when glucose levels seem quite lower than test value
- When hypoglycemia is not well recognised
- When insulin has the biggest effect
- 2 hours after an insulin injection

AST Precautions

- Before using AST, please consult your healthcare professional.
- Do not ignore the symptoms of hyperglycemia or hypoglycemia.
- When the results of the test do not reflect your opinion, retest using the fingertip test. If the fingertip result still does not reflect the way you feel, please consult your healthcare professional.
- Do not rely on the AST results for changing your treatment method.
- The amount of glucose in alternate sites differs from person to person.

Note:

- Results from alternate sites and fingertip samples may differ from each other as there is a time lag for the glucose levels to reach the same value. Use a fingertip for testing if you suffer from hypoglycemia or have experienced hypoglycemic shock or symptoms.
- If the sample drop of blood runs or spreads due to contact with hair or with a line in your palm, do not use that sample. Try puncturing again in a smoother area.

HI and Lo Messages

HI Message

The meter displays results between 1.1–33.3 mmol/L. 'HI' appears when the blood glucose level is greater than 33.3 mmol/L and indicates severe hyperglycemia (much higher than normal glucose levels).

If 'HI' is displayed again upon retesting, please contact your healthcare professional immediately.



Lo Message

'Lo' appears when a test result is less than 1.1 mmol/L and indicates severe hypoglycemia (very low glucose levels).

If 'Lo' is displayed again upon retesting, please contact your healthcare professional immediately.



Target Blood Glucose Ranges

Reminders

Time of day

*Your target ranges
from your healthcare professional*

Before breakfast

Before lunch or dinner

1 hour after meals

2 hours after meals

Between 2 a.m. and 4 a.m.

Expected Values : Normal blood glucose levels for an adult without diabetes are below 5.5 mmol/L before meals and fasting* and are less than 7.8 mmol/L two hours after meals. *Fasting is defined as no caloric intake for at least eight hours.

Reference

American Diabetes Association (Standards of Medical Care in Diabetes – 2021. Diabetes Care), January 2021, vol. 44 (Supplement 1): S15-S33.

Transferring Test Results Using Cable

Test results stored in CareSens N Premier meter can be transferred from the meter to a computer using SmartLog™ software and cable. The meter screen displays 'Pc' when it is connected to the computer using the data cable. For more information, contact your authorised i-SENS sales representative or visit www.i-sens.com.



Note:

Test results can also be transferred wirelessly using Bluetooth. To pair the meter and your smartphone, see page 19.

Meter Memory

The CareSens N Premier meter can save up to 1,000 glucose test results with time and date. If the memory is full, the oldest test result will be deleted and the latest test result will be stored. The meter calculates and displays the averages of total test results, Pre-meal (🍏) test results, Post-meal test (🍷), and Fasting test results (🕒) from the last 1, 7, 14, 30 and 90 days.

Viewing Averages Stored in Memory

Step 1


Press the ◀, ▶ or S button to turn the meter on. The current date and time will be displayed at the bottom of the screen followed by the 1 day average value and the number of the test results saved within the current day.



The number of tests within the current day


Step 2 Viewing Averages

Press the ◀ button to view 7, 14, 30 and 90-day average values and the number of tests performed for the last test period.

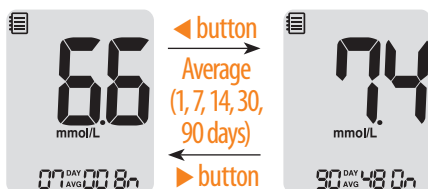


Step 3 Viewing Pre-meal Averages

Repeatedly press the ◀ button to view 1, 7, 14, 30 and 90-day average values and the number of tests performed pre-meals with the 🍏 symbol for the last test period.

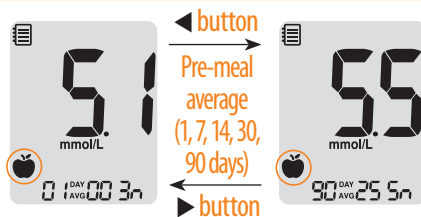


Step 4 Viewing Post-meal Averages

Press the ◀ button to view 1, 7, 14, 30 and 90-day average values and the number of tests performed post-meals with the 🍷 symbol for the last test period.

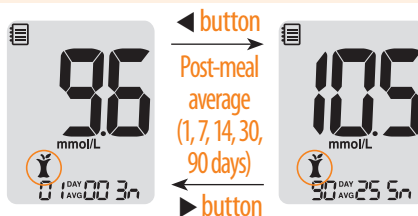

Step 5 Viewing Fasting Averages

Press the ◀ button to view 1, 7, 14, 30 and 90-day average values and the number of tests performed during fasting with the 🚫 symbol for the last test period.

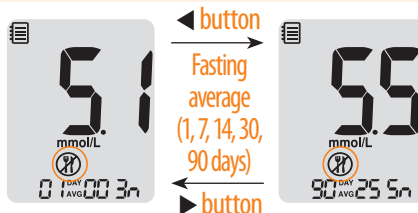


Step 6

Use the ► button to scroll back through the averages seen previously.

Press the S button to turn off the meter.

Note:

The control solution test results saved with the  symbol are not included in the averages.

Viewing Test Results Stored in Memory

Step 1

Press the ◀, ▶ or **S** button to turn the meter on. The current date and time will be displayed on the bottom of the screen followed by the 1 day average value and the number of the test results saved within the current day.




The number of tests within the current day

Step 2



Use the ▶ button to scroll through the test results, starting from the most recent and ending with the oldest.

Press the ◀ button to return to the result seen previously.

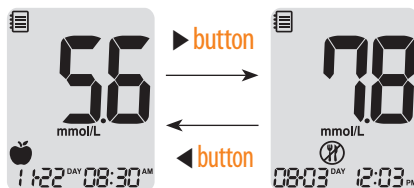After checking the stored test results hold the **S** button to turn off the meter.



Note:

The control solution test results saved with  symbol will be displayed with  symbol when you review the stored test results.

Setting the Alarm Function


Four types of alarms can be set in the CareSens N Premier Meter: one post-meal alarm (PP2 alarm) and three time set alarms (alarm1–3). The PP2 alarm goes off 2 hours after setting the alarm. The alarms ring for 15 seconds and can be silenced by pressing the ◀, ▶ or S button or by inserting a test strip.

Setting the Post-meal Alarm (PP2 alarm)

Step 1 Turning the PP2 alarm On

Without inserting a test strip, press and hold the ◀ button for 3 seconds to set the post-meal alarm.

'PP2', bell (🔔) symbol and 'On' will be displayed. The screen will then automatically change to the memory recall mode. At this time, bell (🔔) symbol, indicating that the PP2 alarm has been set, will be displayed on the screen.

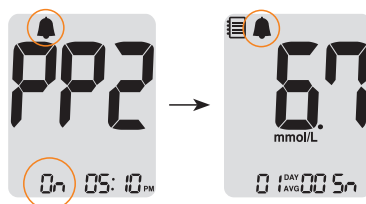


Note:

The PP2 alarm will automatically turn off if the meter's time setting is adjusted to more than two hours before or just past the currently activated PP2 alarm time.

Step 2 Turning the PP2 alarm OFF

To turn off the PP2 alarm, press and hold the ◀ button for 3 seconds.

'PP2', bell (🔔) symbol and 'OFF' will appear on the screen. Then the screen will change automatically to the memory recall mode without bell (🔔) symbol displayed.

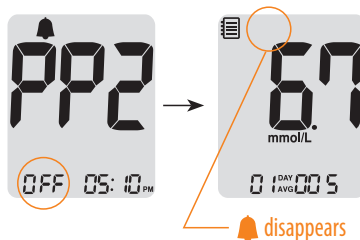

Setting the Time Alarms (alarm 1–3)

Step 1

Without inserting a test strip, press the ◀ and **S** buttons simultaneously for 3 seconds to enter the time alarm setting 'alarm 1' will be displayed while 'OFF' blinks on the screen.



Step 2

On pressing the ▶ button, 'alarm 1' is set and 'On' is displayed on the screen. Press the ▶ button again to cancel 'alarm 1'. 'OFF' will blink on the screen.



Step 3

Press the ◀ button to adjust the time of 'alarm 1'.
A number representing the hour will blink on the screen. Press the ▶ button to set the hour.



Step 4

On pressing the ◀ button, the number indicating the minute will start blinking. Press the ▶ button to set the minute.



Step 5







Press the **S** button to finish and to go to 'alarm 2' setting. Repeat steps 2 to 4 to set the remaining time alarms (alarm 2–3).



Step 6

Press the **S** button for 3 seconds to finish and turn the meter off.

Understanding Error Messages

| | |
|--|---|
|  | <p>A used test strip was inserted. > <i>Repeat the test with a new test strip.</i></p> |
|  | <p>The blood or control solution sample was applied before the  symbol appeared. > <i>Repeat the test with a new test strip and wait until the  symbol appears before applying the blood or control solution sample.</i></p> |
|  | <p>The temperature during the test was above or below the operating range. > <i>Move to an area where the temperature is within the operating range (5–50 °C) and repeat the test after the meter and test strips have reached a temperature within the operating range.</i></p> |
|  | <p>The blood sample has abnormally high viscosity or insufficient volume. > <i>Repeat the test after inserting a new test strip.</i></p> |

Er5

This error message may appear when the wrong blood glucose test strip is used instead of CareSens N blood glucose test strip.
> *Repeat the test with a CareSens N test strip.*

Er6

There is a problem with the meter.
> *Do not use the meter. Contact your authorised i-SENS sales representative.*

Er7

There is a problem with Bluetooth communication
> *Contact your authorised i-SENS sales representative.*

Er8

An electronic error occurred during the test, commonly caused by foreign materials in the meter port (e.g., blood, dirt, dust), and/or physical damage to the meter or strips due to incorrect use or storage.
> *Repeat the test with a new test strip. If the error message persists, contact your authorised i-SENS sales representative.*

Note:

If the error messages persist, contact your authorised i-SENS sales representative.

General Troubleshooting

| Problem | Troubleshooting |
|--|--|
| The display is blank even after inserting a test strip. | <ul style="list-style-type: none">• Check whether the test strip is inserted with the contact bars facing up. Check if the strip has been inserted completely into the test strip port.• Check if the appropriate test strip was used.• Check whether the batteries are inserted with the + side facing up.• Replace the batteries. |
| The test does not start even after applying the blood sample on the strip. | <ul style="list-style-type: none">• Check if the confirmation window is filled completely.• Repeat the test after inserting a new test strip. |
| The test result doesn't match the way you feel. | <ul style="list-style-type: none">• Repeat the test after inserting a new test strip.• Check the expiration date of the test strip.• Perform control solution test. |

Note:

If the problem is not resolved, please contact your authorised i-SENS sales representative.

Performance Characteristics

System Accuracy and Measurement Precision

The performance of CareSens N Premier Blood Glucose Monitoring System has been evaluated in laboratory and in clinical tests.

Accuracy: The accuracy of the CareSens N Premier Blood Glucose Monitoring System (Model GM01AAB) was assessed by comparing blood glucose results obtained by patients with those obtained using a YSI Model 2300 Glucose Analyzer, a laboratory instrument.

The following results were obtained by diabetic patients at clinic centers.

| | |
|-----------------------------|-----------------|
| Slope | 1.0041 |
| Y-intercept | 0.1 mmol/L |
| Correlation coefficient (r) | 0.9933 |
| Number of Samples | 600 |
| Range tested | 1.6–28.4 mmol/L |

System accuracy results for glucose concentration < 5.55 mmol/L

| Within ± 0.28 mmol/L | Within ± 0.56 mmol/L | Within ± 0.83 mmol/L |
|--------------------------|--------------------------|--------------------------|
| 150/186 (80.6 %) | 172/186 (92.5 %) | 186/186 (100.0 %) |

System accuracy results for glucose concentration ≥ 5.55 mmol/L

| Within ± 5 % | Within ± 10 % | Within ± 15 % |
|------------------|-------------------|-------------------|
| 333/414 (80.4 %) | 404/414 (97.6 %) | 414/414 (100.0 %) |

System accuracy results for glucose concentrations between 1.6 mmol/L and 28.4 mmol/L

| |
|--|
| Within ± 0.83 mmol/L and Within ± 15 % |
| 600/600 (100.0 %) |

Precision: The precision studies were performed in a laboratory using CareSens N Premier Blood Glucose Monitoring Systems.

| <i>Within Run Precision</i> | | |
|-----------------------------|-------------|-----------------|
| Blood average | 2.3 mmol/L | SD = 0.1 mmol/L |
| | 4.0 mmol/L | SD = 0.1 mmol/L |
| | 7.1 mmol/L | CV = 3.9 % |
| | 10.4 mmol/L | CV = 3.1 % |
| | 16.2 mmol/L | CV = 2.9 % |

| <i>Between Run Precision</i> | | |
|------------------------------|-------------|-----------------|
| Control average | 2.1 mmol/L | SD = 0.1 mmol/L |
| | 6.7 mmol/L | CV = 3.4 % |
| | 19.7 mmol/L | CV = 4.4 % |

This study shows that there could be variation of up to 4.4 %

Influence Quantities

Packed Cell Volume (Hematocrit)

Packed cell volume evaluation was conducted in various hematocrit levels. The range of hematocrit levels within the acceptance criteria is 15–65 %.

Interferences

The effect of various interfering substances was evaluated in whole blood samples. The presence of the following substances within the given concentrations does not affect blood glucose measurements. Higher concentrations of the substances shown below may cause inaccurate blood glucose results.

| No. | Interferent | Concentration |
|-----|-----------------------------|---------------|
| 1 | Acetaminophen (paracetamol) | 20 mg/dL |
| 2 | Ascorbic acid | 3 mg/dL |
| 3 | Bilirubin | 20 mg/dL |
| 4 | Ceftriazone | 81 mg/dL |
| 5 | Cholesterol | 500 mg/dL |
| 6 | Creatinine | 3 mg/dL |
| 7 | Dopamine | 13 mg/dL |
| 8 | EDTA | 180 mg/dL |
| 9 | Galactose | 60 mg/dL |
| 10 | Gentisic acid | 50 mg/dL |
| 11 | Glutathione (Red) | 17 mg/dL |
| 12 | Hemoglobin | 500 mg/dL |
| 13 | Heparin | 8000 U/dL |
| 14 | Hydrocortisone | 1 mg/dL |
| 15 | Ibuprofen | 40 mg/dL |
| 16 | Icodextrin | 2 mg/dL |

| No. | Interferent | Concentration |
|-----|--|---------------|
| 17 | L-Dopa (L-3,4-dihydroxyphenylalanine) | 5 mg/dL |
| 18 | Maltose | 1000 mg/dL |
| 19 | Mannitol | 600 mg/dL |
| 20 | Methyl-DOPA | 1.5 mg/dL |
| 21 | Pralidoxime Iodide (PAM) | 25 mg/dL |
| 22 | Salicylate | 60 mg/dL |
| 23 | Tolbutamide | 100 mg/dL |
| 24 | Tolazamide | 100 mg/dL |
| 25 | Triglycerides | 3000 mg/dL |
| 26 | Uric acid | 20 mg/dL |
| 27 | Xylose | 25 mg/dL |

User Performance Evaluation

A study evaluating glucose values from fingertip capillary blood samples obtained by 100 lay persons showed the following results:

96.7 % within ± 0.83 mmol/L of the medical laboratory values at glucose concentrations below 5.55 mmol/L, and 99 % within ± 15 % of the medical laboratory values at glucose concentrations at or above 5.55 mmol/L.

Warranty Information

Manufacturer's Warranty

i-SENS, Inc. warrants that the CareSens N Premier Meter shall be free of defects in material and workmanship in normal use for a period of five (5) years. The meter must have been subjected to normal use. The warranty does not cover improper handling, tampering, use, or service of the meter. Any claim must be made within the warranty period.

i-SENS will, at its discretion, repair or replace a defective meter or meter part that is covered by this warranty. As a matter of warranty policy, i-SENS will not reimburse the consumer's purchase price.

Obtaining Warranty Service

To obtain warranty service, you must return the defective meter or meter part along with proof of purchase to your nearest i-SENS sales or customer service representative.

For Australian Consumers

Our goods come with guarantees that cannot be excluded under the Australian Consumer Law. You are entitled to a replacement or refund for a major failure and compensation for any other reasonable foreseeable loss or damage. You are also entitled to have the goods repaired or replaced if the goods fail to be of acceptable quality and the failure does not amount to a major failure.

NZ Distributor:

Pharmaco (NZ) Ltd
Auckland
Phone: 0800 GLUCOSE (0800 458267)
Website:
www.pharmacodiabetes.co.nz

Australian Distributor:

Pharmaco (Australia) Ltd
Chatswood, NSW 2067
Phone: 1800 114 610
Website:
www.pharmacodiabetes.com.au

MEMO

MEMO

MEMO

MEMO

MEMO